

zumax



**SEE BETTER
PERFORM
BETTER**
ZUMAX® Microscopes


enova®

REVOLUTIONIZE YOUR PRACTICE

WITH A ZUMAX® MICROSCOPE

Magnification using the latest in microscope technology gives professionals like you the opportunity to offer a higher standard of care in dentistry. Enova Illumination is proud to be the exclusive U.S. distributor of Zumax® microscopes, bringing you the highest quality and latest innovations available in the USA.



ENHANCED MAGNIFICATION

Increased precision in all aspects of patient care - from diagnosis and treatment, to efficiency and patient outcomes.



INTEGRATED 4K CAMERA

Connect to a monitor or your network for the ability to stream and save high definition photos and videos. Edit videos immediately for training and sharing with consulting physicians. Display real-time imaging for increased patient education.



SUPERIOR ERGONOMICS

Reduce neck and back injuries that can limit a career by not relying on loupes and insufficient illumination. Improve your posture by utilizing heads-up magnification with a monitor. 3D visualization is optional.



DOCUMENTATION AND TRAINING

Microscopes allow you to easily show, discuss and document treatments and procedures, both with patients and other specialists.



HIGH QUALITY OPTICS

On all Zumax® microscopes, all lenses are multi-coated, anti-reflective, apochromatic, and made from hand-selected German glass.



RETURN ON INVESTMENT

Elevate your practice with Zumax® Microscopes, empowering you to confidently perform complex procedures with enhanced visual precision!

Zumax® Dental Microscopes

3200 R2/PRO

FEATURES

- **Joystick-Controlled Continuous Zoom:**
Magnification ranging from ~1-20x
AutoFocus is standard
- **Electromagnetic Clutch:**
One-click magnetic clutch for ease of positioning and smooth movement
- **Integrated 4K Camera:** Seamlessly save, record and stream images/videos to your network and monitor(s)
*3D 4K camera optional
- **Integrated LED Illumination:**
≥105,000 lux at 250mm working distance, multiple spot sizes and light and laser filters included. LED lifetime: 60,000 hours
- **Integrated VarioScope:**
Working range 200mm to 450mm



**Motorized Continuous
Zoom System**



**Electromagnetic
Clutch**



**Integrated 4K
Camera**

Zumax® Dental Microscopes

2380 R2

FEATURES

- **Continuous Zoom Magnification:** Ranging from ~1-20x
- **FreeFlow balancing arm:** Free flowing, ultra smooth movement to simplify positioning
- **Integrated 4K Camera:** Seamlessly save, record and stream images/videos to your network and monitor(s)
*3D 4K camera optional
- **Integrated LED Illumination:**
≥100,000 lux at 250mm working distance, multiple spot sizes and light and laser filters included. LED lifetime: 60,000 hours
- **Integrated VarioScope:**
Variable working distance 200mm-450mm



**Continuous
Zoom System**



**VarioScope Fine
Focus**



**Integrated 4K
Camera**

Zumax® Dental Microscopes

2050

FEATURES

- **Continuous Zoom Magnification:**
Ranges from ~1-20x
- **ErgoPendulum Arm:**
Binoculars remain stationary as microscope pivots right and left
- **Integrated LED Illumination:**
≥75,000 LUX at 250mm working distance, multiple spot sizes and light and laser filters included. LED lifetime: 60,000 hours
- **4K Camera Module (Optional):**
Seamlessly save and display images and video
- **Tabletop Mount:**
Available option



**Continuous
Zoom System**



Integrated LED

Zumax® Dental Microscopes

2350

FEATURES

- **6-Step Magnification:**
Ranging from ~1-24x
- **Integrated LED Illumination:**
≥65,000 LUX at 250mm working distance, multiple spot sizes and light and laser filters included. LED lifetime: 60,000 hours
- **Integrated HD Camera (Optional):** Seamlessly connects to your network
- **Integrated VarioScope:**
Variable working distance 190mm-300mm



Integrated LED



**VarioScope Fine
Focus**

ACCESSORIES

BROKEN INSTRUMENT REMOVAL KIT

Instrument separation is one of the most common complications in the root canal treatment process. The innovative design of the Zumax® Broken Instrument Removal Kit makes instrument removal quick and easy.


Kit Includes:

Extractor Handle and Tip

- 3 Micro-Trephine Burs
- 1 Extractor Handle
- 3 Extractors
- 1 Foreign Forceps



Extractor Tip Options

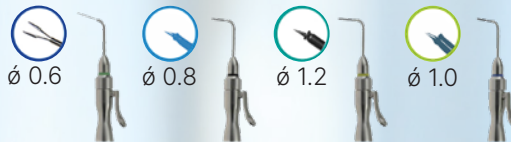


MICRO DENTAL TWEEZER

The Zumax® Micro Dental Tweezer is designed for accurate foreign material removal - no matter the angle or complexity.

Kit Includes:

- 1 360 Tweezer Head
- 1 Handle
- 1 Cleaning Adapter

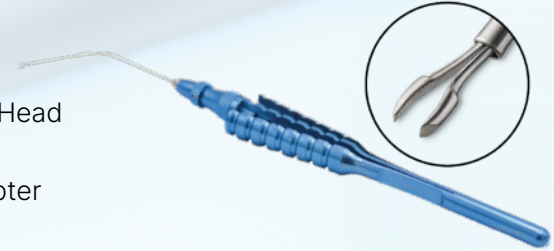

MOBILE ADAPTERS

Connect mobile devices to Zumax® microscopes to allow you to instantly share your images through live-streaming, text messages, or email.



Contact your Enova sales rep today to learn more about Zumax microscopes and accessories.

Visit www.goenova.com, call 651.236.8857, or email sales@goenova.com.

ENOVA SATISFACTION GUARANTEE

Use our products and love them. If not, we'll take them back. It's that simple.
For full warranty specifications visit www.goenova.com/product-support

Connect with us!



To learn more about how Enova can help you see better, so you can perform better,
visit **www.goenova.com**, call 651.236.8857, or email sales@goenova.com

World Headquarters

Enova Illumination
6305 Sandburg Road, STE 100
Minneapolis, MN 55427

ISO 13485: 2016 Certified QMS
MDR Class 1 Certified

F20578B

© 2023 Enova Illumination, Inc. All trademarks are the property of Enova Illumination, Inc. unless otherwise noted.
Oakley® and Flak® are trademarks of Oakley, Inc.
Zumax® is a registered trademark of Zumax Medical Co., LTD

